Vital Signs™ Gas sampling line, disposable

REF

73318-HEL Vital Signs™ Anesthesia Gas Sampling Line, Disposable, 2 m/7 ft
73318-EN Vital Signs™ Anesthesia Gas Sampling Line, Disposable, 2 m/7 ft
73318-FR Vital Signs™ Anesthesia Gas Sampling Line, Disposable, 3 m/10 ft
73319-DE Vital Signs™ Anesthesia Gas Sampling Line, Disposable, 3 m/10 ft

REF	Catalogue number	MD	Medical device
LOT	Lot number		Consult instructions for use
QTY	Quantity		Follow instructions for use
	Expiration date		Unique device identifier
	Manufacturer		Single use
	Date of manufacture		Not made with natural rubber latex.
	Storage Temperature		Indicates the authorised representative in Switzerland
	The UK Conformity Assessed marking declares compliance with the applicable UK legislation.		Cautlon: U.S. Federal law restricts this device to sale by or on the order of a licensed healthcare practitioner
	The CE marking declares compliance with the applicable EU legislation.		


EN English

INSTRUCTIONS FOR USE

INTENDED USE:
The gas sampling lines are used with GE HealthCare Compact Airway modules and CARESCAPE™ Respiratory modules for monitoring CO₂, O₂, N₂O and anesthetic agents and with E-minic® for monitoring CO₂.
These devices are indicated for use by healthcare professionals only.

INDICATIONS FOR USE:
Gas Sampling Lines are indicated to sample exhaled gas in intubated or non-intubated patients requiring expired gas monitoring.

ENVIRONMENT FOR USE:
Hospital, sub-acute, surgical and centers.

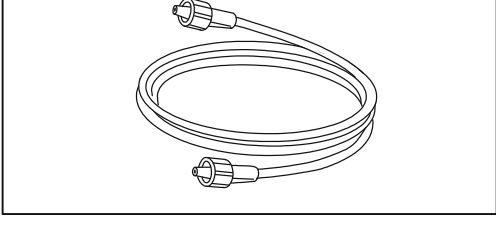PATIENT RISK GROUP:
Pediatric, Adult.

NECESSARY USER QUALIFICATIONS:
The end-user should receive training on the proper use of the device and be familiar with the indications for use, warnings, and cautions as stated in the IFU.
Rx Only: Federal law (USA) restricts this device to sale by or on the order of a licensed healthcare practitioner.

EXPECTED CLINICAL BENEFITS:
Capnography provides early detection and warning of changes in a patient's ventilatory status. This enables to clinicians to adjust ventilation dosage per applied therapeutic intervention to prevent hypoventilation. EICO₂ monitoring provides information on endotracheal tube placement, the effectiveness of chest compressions and the return of circulation.
Accurate monitoring of anesthesia gases ensure the patient receives the correct dose to maintain appropriate levels of anesthesia and prevent intraoperative awareness. Monitoring the concentration of inspired and expired anesthetic levels, clinicians can titrate the dose precisely, preventing the potential for an anesthetic overdose.

CONTRAINDICATIONS:
None.

WARNINGS:
- Do not place anything on sampling line that can obstruct flow.
- Check sample line regularly during use to ensure that no kinks are present. Kinked tubing may cause loss of or inaccurate readings.
- Check connections. Loose connections may compromise ventilation or cause inaccurate gas measurements.
- Do not sterilize or clean device as it may affect the performance of the device and cause damage to the monitor.
- When used with a closed suction catheter (CSC) and/or heat moisture exchanger (HME), place the Sampling Line and Airway Adapter between the circuit wye and the other devices.

MATERIAL:
TPE (Thermoplastic elastomer)

RESIDUAL RISKS:
Refer to warnings and cautions as applicable.

CAUTIONS:
- Do not place anything on sampling line that can obstruct flow.
- Check sample line regularly during use to ensure that no kinks are present. Kinked tubing may cause loss of or inaccurate readings.
- Check connections. Loose connections may compromise ventilation or cause inaccurate gas measurements.
- Do not sterilize or clean device as it may affect the performance of the device and cause damage to the monitor.
- When used with a closed suction catheter (CSC) and/or heat moisture exchanger (HME), place the Sampling Line and Airway Adapter between the circuit wye and the other devices.

MATERIAL:
TPE (Thermoplastic elastomer)

RESIDUAL RISKS:
Refer to warnings and cautions as applicable.

INSTRUCTIONS FOR USE:
1. Inspect packaging and product prior to use. Do not use if the integrity of the packaging has been compromised or product is damaged.
2. Connect one end of the sampling line to the patient circuit and the other end to the gas input port on the anesthesia gas monitor.
3. Twist the connector clockwise until you feel resistance.
4. Check all connections for leaks according to standard clinical procedures.
5. Ensure monitor gas valves and EICO₂ waveforms are displayed on monitor.
6. Do not reuse sampling line. Dispose of sampling line if it becomes dirty or damaged.

STORAGE:
Storage temperature -20°C (-4°F) to +50°C (+122°F).

SAFE DISPOSAL:
Decontaminate and dispose of all potentially biohazardous material. Dispose of device in accordance with local, state or national regulations.

INCIDENT REPORTING:
Notice to the user that any serious incident that has occurred in relation to the device should be reported to the manufacturer at productquality@myairlife.com and to the competent authority of the Member State in which the user and/or patient is established.

TECHNICAL DESCRIPTION:
For technical description, see the host device user guide or manual, or contact your local sales representative. Not made with natural rubber latex.

ORDERING INFORMATION:
For further information, please contact your local sales representative or visit www.myAirLife.com.

	AirLife Finland Oy Kuortaneentie 2 FI-02015 Helsinki, Finland www.myAirLife.com
	2797
	0086

UK Responsible Person:
MT Promed Consulting Ltd.
First Floor, Park Central
1201 Park East Street
Oxford OX1 1UD
United Kingdom

Australian Sponsor:
AirLife Australia Holdings Pty Ltd.
PO BOX 97
NORTH RYDE NSW 1670
Australia

CH	REP
	MedEnvoy Switzerland Göthelstrasse 28 8302 Zug Switzerland www.medivital.com

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