(en)

1101 Lucas Ave. • Suite 200
St. Louis, MO 63101-1159 USA
Phone: 314-781-5700 • Fax: 314-644-4262
email: airQ@cookgas.com • www.cookgas.com

air-Q SP® 3/3G
INTUBATING LARYNGEAL AIRWAY
Instructions for Use

The air-Q SP 3/3G (Self-Pressurizing) Airway is an enhanced version of the standard air-Q masked laryngeal airway. As such, it is indicated as a primary airway device when an oral endotracheal tube (OETT) is not necessary or as an aid to intubation in difficult airway situations.

The air-Q SP 3/3G represents a revolution in airway device design. It functions similarly to the original air-Q in that it retains the soft perimeter mask cuff. This gives the mask the ability to change its size and shape depending on each patient's pharyngeal anatomy, improving fit. Like the original, it can also be used as a conduit for intubation in difficult airway situations. It differs in that it does not contain an inflation valve and so does not inflate normally. Instead, the air-Q SP 3/3G incorporates a large aperture between the internal volumes of the breathing tube and the peripheral cuff. This aperture allows fluid communication between the tube and cuff to "self-pressurize", causing inflation of the cuff during positive pressure ventilation.

This product is to be used by trained personnel only.

Available in Single Use Only

Recommendations:

Size	IBW	Max. OETT	Mouth Opening ¹	← → ²	Max OG Tube ³
5	>80-kg	9.0 mm	25 mm	20 cm	18
4	60-80 kg	8.0 mm	23 mm	18 cm	18
3	30-60 kg	7.0 mm	20 mm	16 cm	16
2	17-30 kg	5.5 mm	17 mm	14 cm	12
1.5	7-17 kg	5.0 mm	14 mm	11 cm	10
1.0	4-7 kg	4.5 mm	11 mm	9 cm	9
0.5	2-4 kg	4.0 mm	8 mm	7 cm	8
0	≤ 2 kg	3.0 mm	5 cm	6 cm	7

- Minimum mouth opening for insertion.
- Distance from the external edge of the airway tube to the internal ventilatory opening.
- Maximum OG size (airQ SP 3G only).

- air-Q SP 3/3G Placement Procedure**
- The procedure below is intended as a guide. Many techniques can be successfully used to properly place the air-Q SP 3/3G into the pharynx.
- Lubricate the external surface including the mask cavity ridges.
 - Open the patient's mouth and elevate the tongue. Elevating the tongue lifts the epiglottis off the posterior pharyngeal wall and allows the air-Q SP 3/3G easy passage into the pharynx. A mandibular lift is especially recommended. A tongue blade placed at the base of the tongue also works well for this purpose.
 - Place the front portion of the air-Q SP 3/3G mask between the base of the tongue and the soft palate at a slight forward angle, if possible.
 - Pass the air-Q SP 3/3G into position within the pharynx by gently applying inward and downward pressure, using the curvature of the air-Q SP 3/3G mask and airway tube as a guide. Simply rotate the air-Q SP 3/3G forward and inward. Minimal manipulation may be necessary to turn the corner into the upper pharynx. Continue to advance until fixed resistance to forward movement is felt. Correct placement is determined by this resistance to further advancement. Some users place the back of the left index finger behind the mask, flexing the finger forward to help guide the mask around the corner into the pharynx. Once the mask has negotiated the turn, the left hand is then used to do a mandibular lift while exerting downward and inward pressure on the air-Q SP 3/3G with the right hand during final advancement into the pharynx. This technique seems to be easy to learn and is particularly successful.
 - Tape the air-Q SP 3/3G in place.
 - Check the air-Q SP 3/3G connector to ensure it is fully engaged within the airway tube, and attach the connector to the appropriate breathing device. Check for adequate ventilation.
 - Place a bite block between the patient's teeth. Keep the bite block in place until the air-Q SP 3/3G is removed.

air-Q SP 3/3G Intubation Procedure

The air-Q SP 3/3G by Cookgas® LLC is intended not only to be an outstanding airway for general use, but also to be a simple and reliable tool for intubation of the trachea with OETTs. Due to its patented design, standard OETTs (sizes 9.0 mm - 3.0 mm) can be easily passed through the air-Q SP 3/3G and into the trachea. Further, the air-Q SP 3/3G can be easily removed following intubation with the aid of the patented air-Q Removal Stylet, also by Cookgas® LLC. The following procedure for intubation is only intended as a guide. Many techniques can be successfully used for tracheal intubation using the air-Q SP 3/3G.

- Prior to intubation, the laryngeal musculature and vocal cords must be relaxed, either by an aerosolized local anesthetic or with the aid of a muscle relaxant.
- Pre-oxygenate.
- Prepare the appropriately-sized OETT by completely deflating the OETT cuff and lubricate well. It is important to deflate the OETT cuff completely to allow the OETT to slide easily within the air-Q SP 3/3G.

- Disconnect the air-Q SP 3/3G from the ventilation device and remove the air-Q SP 3/3G connector. This can be easily done by squeezing the air-Q SP 3/3G tube between the index finger and thumb just distal to the connector with one hand, then rocking the air-Q SP 3/3G connector back and forth while pulling the connector outward away from the airway tube with the other.
- Insert the previously deflated and lubricated OETT through the air-Q SP 3/3G to a depth of approximately 6 - 20 cm, depending on the air-Q SP 3/3G size. This will place the distal tip of the OETT at or just proximal to the opening of the air-Q SP 3/3G airway tube within the mask cavity. It is very important to lubricate the OETT and the air-Q SP 3/3G airway tube completely to ensure easy passage of the OETT through the air-Q SP 3/3G.
- The following suggestions for advancement of the OETT are intended as a guide. Many techniques can be successfully used to further advance the OETT into the trachea.

- CAUTION:** Always check for adequate ventilation and oxygenation following placement of the OETT.
- Fiber Optic Technique:** Using a Fiber Optic Endoscope, pass the scope through the OETT and into the trachea under direct visualization. Stabilize the Fiber Optic Endoscope and pass the OETT through the laryngeal inlet and into the proximal trachea, using the scope as a guide. Check the position of the OETT with direct visualization of the tracheal carina. Remove the Fiber Optic Endoscope. Add a small amount of air to the OETT cuff, then replace the OETT connector if necessary. Check for adequate ventilation. (If epiglottic intrusion or down-folding is seen during visualization, the air-Q SP 3/3G usually need not be completely removed.) Perform an external mandibular lift to elevate the epiglottis and pass the endoscope beneath the epiglottis and into the trachea followed by the OETT.
 - Stylet Technique:** Using an appropriate coude-tipped intubating stylet or a lighted stylet, pass the intubation stylet through the OETT within the air-Q SP 3/3G, through the laryngeal inlet and into the trachea. Pass coude-tipped stylets with tip pointing upward (anterior). By gently placing the fingers of the left hand over the cricoid area of the patient's throat, the stylet can usually be felt as a scraping or rubbing sensation as it passes through the cricoid ring. If properly positioned, the lighted stylet will also produce a bright yellow/red illumination over the cricoid area. Once the stylet passes into the trachea, simply advance the OETT over the stylet, through the laryngeal inlet and into the trachea, using the intubation stylet as a guide. Add a small amount of air to the OETT cuff, replace the OETT connector and check for adequate ventilation.

NOTE: If the OETT fails to advance over the stylet into the trachea, it is usually helpful to rotate the OETT counterclockwise while passing the OETT. If this fails, try again with a smaller size OETT.

air-Q SP 3/3G Removal Procedure

Removing the air-Q SP 3/3G following OETT intubation is easily accomplished with the aid of the air-Q Removal Stylet by Cookgas® LLC. The air-Q Removal Stylet consists of an adapter connected to a rod. The adapter is tapered from bottom to top, and has horizontal ridges and vertical grooves. The taper allows the stylet to accommodate multiple OETT sizes. The ridges engage the OETT in a firm, secure grip, giving the user control of the OETT during the air-Q SP 3/3G removal process. The grooves allow spontaneously breathing patients unimpeded air passage within the OETT during removal of the air-Q SP 3/3G. By immobilizing and exerting an inward stabilizing force on the OETT, the air-Q SP 3/3G Removal Stylet allows for the swift removal of the air-Q SP 3/3G without dislodging the previously-placed OETT from the patient.

- Remove the OETT connector from the OETT.
- Squeeze the proximal portion of the OETT between the index finger and the thumb, leaving enough room for the adapter portion of the stylet to enter the proximal opening of the OETT. Alternatively, squeeze the proximal end of the air-Q SP 3/3G airway tube, trapping the OETT inside.
- Insert the tapered end of the air-Q Removal Stylet into the proximal OETT (the long axis should be in the 12 o'clock - 6 o'clock position) until the adapter fits snugly within the OETT.
- For larger sizes, (2.0 - 4.5), with firm inward pressure, rotate the stylet adapter in a clockwise direction (into the 3 o'clock - 9 o'clock position) until the adapter firmly engages the OETT. For smaller sizes (0 - 1.5) simply push the stylet firmly into the OETT. Please practice this a few times prior to attempting on a patient.
- Deflate and lubricate the pilot balloon on the OETT prior to withdrawing the air-Q SP 3/3G. Reinflate the OETT following air-Q SP 3/3G removal.
- While exerting an inward stabilizing force on the stylet, slowly withdraw the air-Q SP 3/3G outward over the stylet rod.
- For larger sizes (2.0 - 4.5) pass the stylet through and through. For smaller sizes (0 - 1.5) remove the stylet from the proximal end of the air-Q SP 3/3G while stabilizing the OETT at the mouth. Discard single-use air-Q SP 3/3Gs following use.
- Reposition the OETT to the proper depth within the patient, if needed, and then tape into place.
- Replace the OETT connector within the OETT. Inflate the OETT, if needed, and attach to an appropriate breathing device. Check for adequate ventilation.

- Cautions/Warnings**
- Inspect all air-Q SP 3/3G devices prior to use. Discard defective devices.
 - Do not use sharp instruments on or near the air-Q SP 3/3G.
 - Confirm that the air-Q SP 3/3G size matches the connector size prior to use.
 - The air-Q SP 3/3G connector is removable and airway disconnect is possible. Take standard precautions to minimize the chance of disconnect.

- Confirm complete connector engagement within the airway tube prior to use.
- The air-Q SP 3/3G connector is removable and airway disconnect is possible. Take standard precautions to minimize the chance of disconnect.
- Do not use excessive force during air-Q SP 3/3G placement or removal.
- Immediately check for adequate ventilation following placement.
- If airway problems occur and persist, remove the air-Q SP 3/3G and establish an effective airway by another method. Back-up means for ventilation should be readily available.
- air-Q SP 3/3G connectors may dislodge during use following lubrication. Clean the breathing tube and connector thoroughly with alcohol prior to re-use.
- Supralaryngeal Airways, including the air-Q SP 3/3G, do not fully protect the patient from aspiration.
- Recheck airway position and patency following all changes in the patient's head or neck position.
- Supralaryngeal Airways are potentially flammable in the presence of lasers and electrical cautery.
- Placement and maintenance of a bite block is recommended during air-Q SP 3/3G use.
- Re-use of single-use devices may lead to mechanical malfunction and potential microbiological contamination. **Discard all single-use air-Q SP 3/3G's following use. DO NOT REUSE.**
- The single-use air-Q SP 3/3G has been sterilized utilizing Ethylene Oxide, a known carcinogen.

Contraindications

The air-Q SP 3/3G is contraindicated in patients at high risk for regurgitation and/or aspiration. This includes, but is not limited to, patients undergoing major thoracic or abdominal surgery, patients who are non-fasted, morbidly obese, pregnant > 14 weeks, or suffer from delayed gastric emptying or esophageal reflux. Users must weigh the benefits of emergency airway needs with the potential risk of aspiration in these patients. air-Q SP 3/3Gs should be used in unconscious or topically anesthetized patients only.

Adverse Effects

Previously reported adverse events with masked laryngeal airways include: sore throat, aspiration, regurgitation, vomiting, bronchospasm, gagging, hiccup, coughing, transient glottic closure, airway obstruction, laryngeal spasm, retching, breath holding, arytenoid dislocation, trauma and/or abrasion to the epiglottis, larynx, pharynx, uvula, hyoid and tonsils, tongue cyanosis, lingual nerve, vocal cord and hypoglossal nerve paralysis, tongue macroglossia, parotid gland swelling, dry mouth, dysphagia, feeling of fullness, mouth ulcer, dysarthria, dysphonia, hoarseness, stridor, pharyngeal ulcer, pulmonary edema, laryngeal hematoma, head and neck edema, myocardial ischemia and dysrhythmia.

Warranties

Cookgas® LLC agrees to warrant the disposable air-Q SP 3/3G for a period of 30 days following the invoice date. Warranty covers materials and manufacturing defects provided that the airway is used according to the procedures outlined in the Instructions For Use (IFU) manual. Warranty is valid only following purchase from authorized distributors.

Cookgas® LLC disclaims all other warranties whether expressed or implied.

Patents U.S.
5,937,860
US 6,422,239 B1
US 6,892,731 B2
US 7,331,347 B2
US 7,900,632 B2
US 9,320,864 B2

Patents CAN
US 6,705,321 B2
US 7,357,845 B2
US 7,780,900 B2
US 7,784,464 B2
US 7,934,502 B2
US 8,622,060 B2
US 8,631,796 B2
US 10,729,866 B2

Patents UK
GB2324040B
GB2407293B
GB2405589B
GB2357437B

Patents JAPAN
5296058

Other USA & Foreign Patents Pending

 **ATTENTION**
SEE DIRECTIONS FOR USE

 Latex Free

 Keep Dry

 Keep Away From Sunlight

 Rx Only
By Prescription Only

 **EXP**
Expiration Date

 Single Use
DO NOT REUSE

 MR Safe

 Do Not Sterilize

STERILE

EO

Distributed exclusively by:

For Ordering Information, Contact:

2710 Northridge Dr. NW, Ste A
Grand Rapids, MI 49544 USA
PH. 800-433-2797
www.Sun-Med.com

EC



REP

M T Promedt Consulting
Altenhofstrasse 80
D-66386 St. Ingbert • Germany
Tel. +49(0) 6894.581020 • Fax +49(0).6894.581021
email: info@mt-procons.com
www.mt-procons.com

Recommended insertion Technique

- Anbefalet indføringsteknik
- Aanbevolen inbrengtechniek
- Technique d'insertion recommandée
- Empfohlenes Einführungsverfahren
- Συνιστώμενη τεχνική εισαγωγής
- Tecnica di inserimento consigliata
- Anbefalt innføringsteknikk
- Zalecana technika wsuwania
- Técnica de Inserção Recomendada
- Рекомендованная техника введения
- Técnica de inserción recomendada
- Rekommenderad införselteknik
- Önerilen Sokma Tekniği
- 推荐插入方法

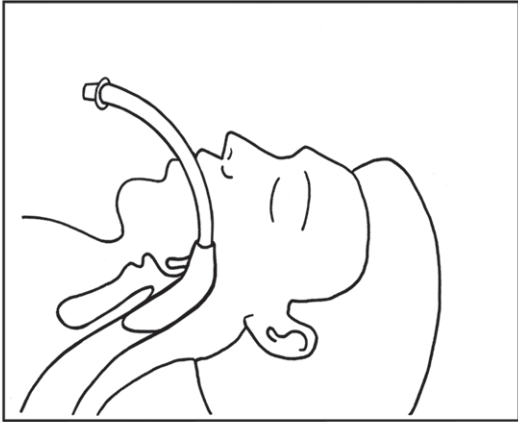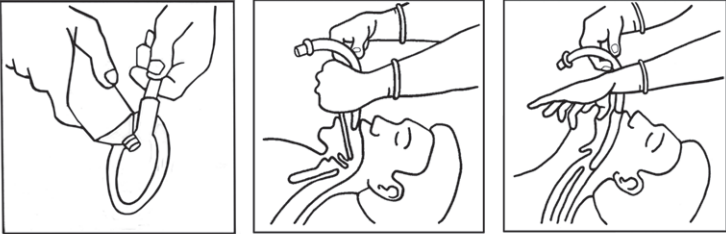Recommended Depth of Insertion Range

- Anbefalet område for indføringsdybde
- Aanbevolen inbrengdieptebereik
- Profondeur recommandée de la portée d'insertion
- Empfohlener Einführtiefenbereich
- Συνιστώμενο βάθος εύρους εισαγωγής
- Gamma di profondità di inserimento raccomandata
- Anbefalt område for innføringsdybde
- Zalecany zakres głębokość wsuwania
- Intervalo de Profundidade de Inserção Recomendada
- Рекомендованный диапазон глубины введения
- Profundidad recomendada del rango de inserción
- Rekommenderat djup för införselspann
- Önerilen Sokma Aralığı Derinliği
- 推荐插入深度的范围

LENGEND

1. Incisors

2. Tongue

3. Hyoid Bone

4. Vocal Cords/folds

5. Thyroid Cartilage

6. Trachea

7. Esophagus

8. Laryngeal Inlet

9. Epiglottis

10. Arytenoid Cartilages

11. Sealing Mechanism

12. Ventilatory Opening

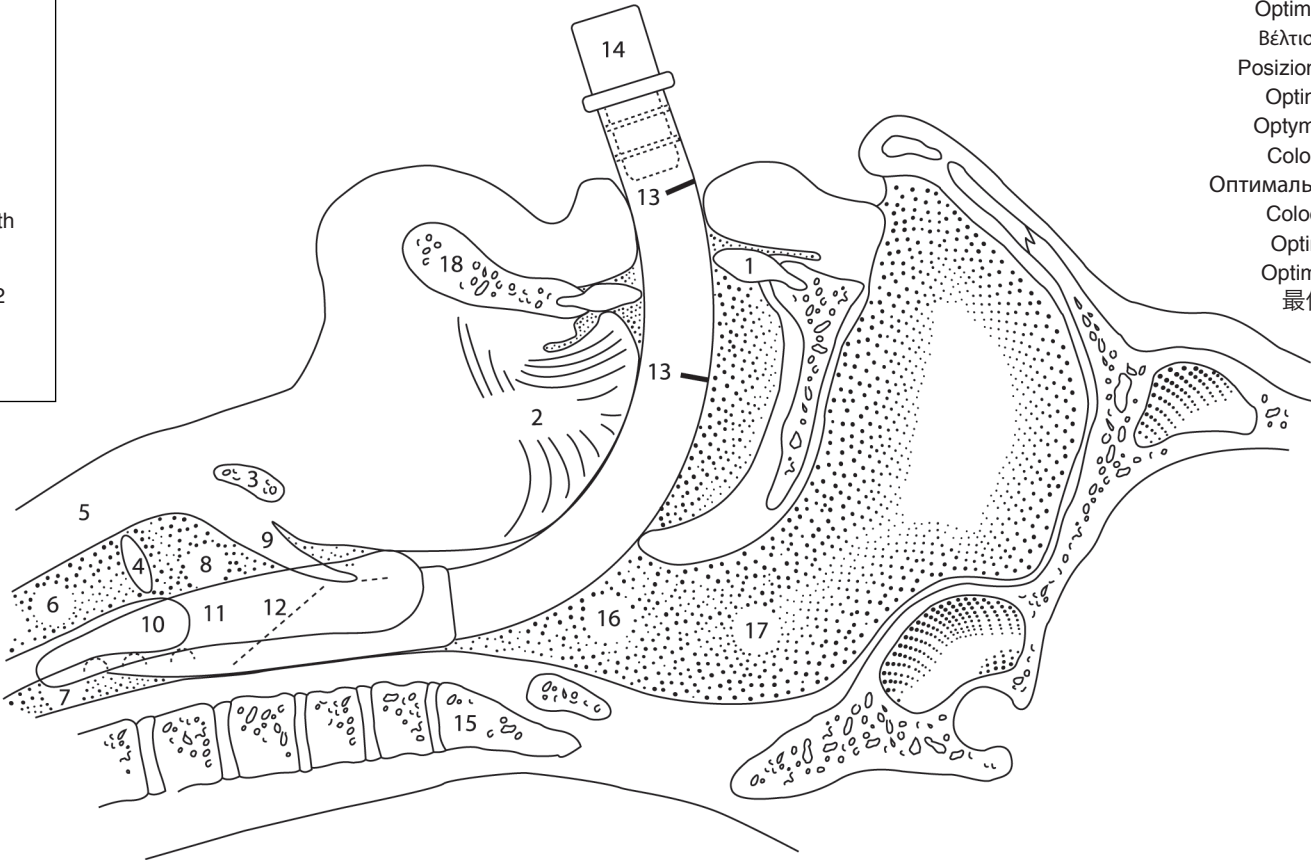
13. Recommended Depth of Insertion Marks

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