


REF 9442**UNITS - 1****Rx Only****CE 0482**

Transport Kit with Adult Elongated Non-Rebreather Mask with Safety Vent

CAUTION: Use of sterilization or cleaning solutions may leave this device nonfunctional.

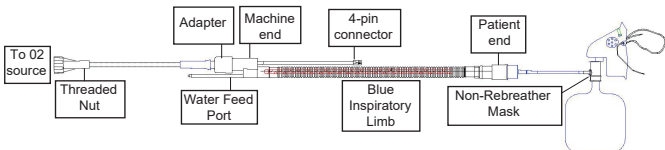

Application: Transport kit is intended to be used in conjunction with Anapod™ Humi-Therm Controller. This kit is for patient transport when using a Non-Rebreather Mask in pediatric and adult applications. Products are intended for use in delivering heated humidification with flow rates up to 12 LPM. All products are to be used in accordance with the applicable prescription.

Patient Transport Setup (Non-Rebreather)

1. Remove Non-Rebreather Mask and 6' O₂ supply tubing from package.
2. Connect 6' O₂ supply tubing to the oxygen source on transport gurney, attach adapter end of oxygen tubing to the inspiratory limb of circuit.
(See Illustration)
3. Attach patient end of inspiratory limb of circuit to Non-Rebreather mask.
(See Illustration)
4. Set oxygen to the desired flow rate and check for flow throughout the device.
5. Place Non-Rebreather Mask on patients face and gently pull both ends of the strap until the mask is secure.
6. Transport patient to the PACU per facility protocol.

PACU Setup

1. Disconnect O₂ supply tubing from oxygen source on transport gurney and connect to wall flow meter in PACU.
2. Connect 4-pin plastic connector on inspiratory limb to Anapod™ controller cable.
3. "Power On" Anapod™ controller and set to desired temperature setting.
4. If there is less than 6" of condensation visible within the inspiratory limb (blue circuit tubing), then add more water. Sterile water is added by attaching a 20ml syringe to the water feed port. If you accidentally overfill the circuit, simply drain excess water at the end of the circuit.

Cautions:

Make sure all fittings are tightly connected before connecting the circuit to the patient. For proper humidification function and maximum safety, make sure the wick is wet as referenced above.

Warning: If breathing circuit is exposed to any form of mechanical compression or damage, replace circuit to ensure localized over heating does not occur. Do not use HME/HMEF with the Anapod™ Humi-Therm circuit.

Packaged clean Nonsterile

Westmed

5580 S. Nogales Highway
Tucson, Arizona 85706 USA
Phone: 800-975-7987
Fax: 520-294-6061
www.westmedinc.com

EC REP

MT Promedt Consulting GmbH
Altenhofstr. 80
66386 St. Ingbert, Germany
Patent Pending
Assembled in Mexico
LABEL PN 75532, Rev. 02



LOT