

REF 001268U

AirLife™

Pediatric Oxygen Mask

Vinyl, Under the Chin Style, 3 in 1 Mask with Safety Vent and 7 foot (2.1 m) Crush Resistant Oxygen Tubing and U/Connect-It Adapter

REF

Catalogue number

R_X ONLY

U.S. Federal law restricts this device to sales by or on the order of a physician

LOT

Lot Number

Clean, Ready to Use

Clean, Ready to Use

QTY

Quantity



Manufacturer



Not made with natural rubber latex



Use-by Date



Contains Phthalates



Do not re-use Single Use




Consult instructions for use



Caution



Temperature Limit

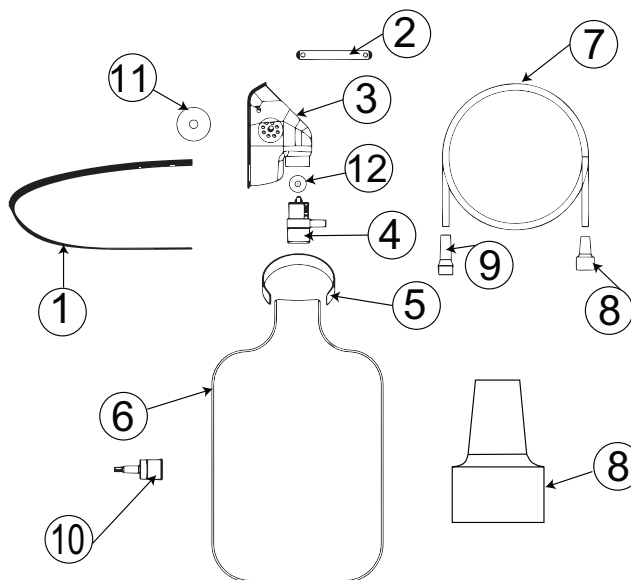

EN

English

Oxygen Mask consists of:

- ① Elastic Strap
- ② Noseclip
- ③ Pediatric Oxygen Mask
- ④ Oxygen Bag Connector
- ⑤ Green Tape
- ⑥ Oxygen Bag
- ⑦ 7 feet (2.1m) Crush Resistant Oxygen Tubing
- ⑧ U/Connect-It Adapter (Threaded Adapter)
- ⑨ Vinyl Tip Adapter
- ⑩ Oxygen Connector
- ⑪ Rebreather Gasket Synthetic
- ⑫ Gasket Synthetic Polyisoprene

AU REF

AirLife Australia Holdings Pty Ltd
PO Box 97
North Ryde BC, NSW, 1670
Australia



AirLife
2710 Northridge Dr. NW, Suite A
Grand Rapids, MI 49544 USA
www.myAirLife.com
Made in Mexico

Intended Use

Intended for oxygen to be delivered to the patient at a prescribed concentration and/or flow rate. Oxygen tubing is connected to an oxygen source on one end and a mask on the other end to deliver oxygen to the patient.

Directions for Use

Mask is assembled as a high concentration, nonbreather mask with safety vent.

1. Attach vinyl tip of oxygen tubing to mask. Attach U/Connect-It Adapter of oxygen tubing directly to threaded oxygen outlet. (Separate tubing adapter may be used but it is not required.) Push and twist to seal. Adjust liter flow as prescribed by physician.
2. Place mask over patient's face covering both nose and mouth.
3. Slip strap over patient's head and adjust for comfort and fit.
4. Pinch nose clip to provide a seal around the nose.

To use as a rebreather, remove flexible disks located on side port and inside of mask on top of reservoir bag. Proceed with steps 1–4.

To use as a medium concentration mask, remove reservoir bag, connector with baffle, and flexible disc on side port. Then, insert oxygen connector (supplied loose in bag) and reconnect tubing. Proceed with steps 1–4.



Caution

When using mask with reservoir bag, make sure bag is inflated before attaching to patient.



Warning

- When used as a nonrebreather mask, assure adequate and uninterrupted oxygen flow to patient. If oxygen flow is inadequate or interrupted, patient may not be able to inhale, or may not be able to inhale adequate volume.
- The oxygen connector in mask should not be connected directly to any endotracheal or tracheostomy tube inlet.
- Avoid kinking or occluding tubing.

Storage

Recommended normal storage temperature 0°C/+32°F to +50°C/+122°F.

Contains Phthalates (DEHP). Risks and Precautionary

Measures related to Phthalates: This instruction pertains to the phthalate symbol marked on the device or its packaging. If this device is used for the treatment of children, or treatment of pregnant or nursing women; please note that the following types of procedures may increase the risk of exposure to phthalates: Exchange transfusion in neonates, total parenteral nutrition in neonates, multiple procedures in sick neonates, haemodialysis in peripuberal males, male foetus and male infant of pregnant women, and lactating women; and massive blood infusion into trauma patients. Although these procedures have the potential for increased risk of exposure, conclusive evidence of human health risks has not been established. As a precautionary measure, to reduce the potential for unnecessary exposures to phthalates, the product must be used in accordance with the instructions for use, and practitioners should refrain from using this product beyond the period of time the product is medically necessary or needed.

Dispose of all materials in accordance with local, state, and federal regulations. Decontaminate and dispose of all potentially biohazardous material.

Any serious incident that has occurred in relation to the device should be reported to the manufacturer at productquality@myAirLife.com.

AirLife®