


CONVENTIONAL MACINTOSH ENGLISH LARYNGOSCOPE BLADES

LATEX FREE | REUSABLE

The AirLife® English profile conventional laryngoscopes are available in both Macintosh and Miller designs. Constructed of 303/304 surgical stainless steel, these blades come with either LED or clear lamps. The LED lamp emits a cool white illumination.

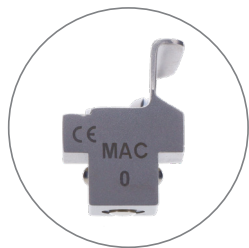


MACINTOSH BLADES WITH LED LAMPS

- Macintosh laryngoscope design is the predominant choice among curved blades
- Flange extends all the way down to distal tip
- Soft matte finish virtually eliminates reflection and glare
- Cool white LED illumination delivers 35,000 hours of use; durable performance plus low current use increases battery life

ENGLISH PROFILE

ITEM	DESCRIPTION	SIZE	A	B	PK
5-5051-00	Neonate	0	80 mm	36.6 mm	1
5-5051-01	Infant	1	95 mm	38.6 mm	1
5-5051-02	Child	2	115 mm	42.3 mm	1
5-5051-03	Medium Adult	3	134 mm	42 mm	1
5-5051-04	Large Adult	4	155 mm	42.3 mm	1
5-5051-05	X-Large Adult	5	175 mm	45.7 mm	1

BLADE MEASUREMENTS

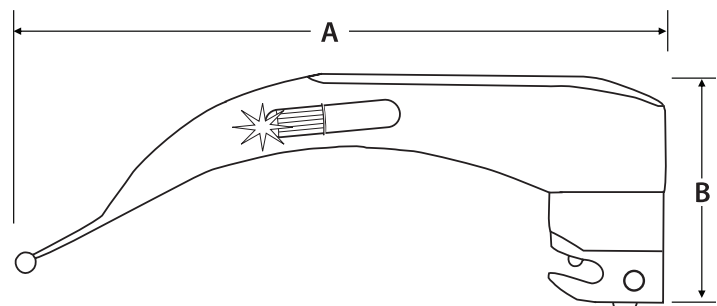

CONVENTIONAL MILLER ENGLISH LARYNGOSCOPE BLADES

LATEX FREE | REUSABLE

The AirLife® English profile conventional laryngoscopes are available in both Macintosh and Miller designs. Constructed of 303/304 surgical stainless steel, these blades come with cool white LED lamps.




CONVENTIONAL MILLER ENGLISH LARYNGOSCOPE BLADES

- Miller laryngoscope design is the predominant choice among straight blades
- Soft matte finish virtually eliminates reflection and glare
- LED and halogen illuminations available
- Cool white LED illumination delivers 35,000 hours of use; durable performance plus low current use increases battery life
- Clear halogen lamp provides bright illumination

CLEAR HALOGEN LAMP LARYNGOSCOPE BLADES

ITEM	DESCRIPTION	SIZE	A	B	PK
5-5061-00	Neonate	0	77 mm	9 mm	1
5-5061-01	Infant	1	102 mm	9 mm	1
5-5061-02	Child	2	153 mm	11 mm	1
5-5061-03	Medium Adult	3	195 mm	11 mm	1
5-5061-04	Large Adult	4	205 mm	16 mm	1
5-0235-01	Replacement Lamp	0, 1	—	—	10
5-0235-04	Replacement Lamp	2, 3, 4	—	—	10

BLADE MEASUREMENTS

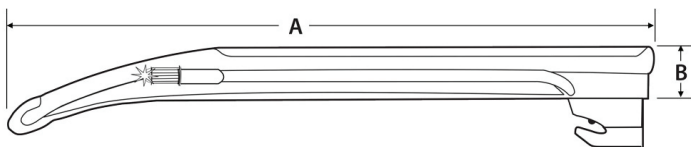

LED LAMP LARYNGOSCOPE BLADES

ITEM	DESCRIPTION	SIZE	A	B	PK
5-5861-00	Neonate	0	77 mm	9 mm	1
5-5861-01	Infant	1	102 mm	9 mm	1
5-5861-02	Child	2	153 mm	11 mm	1
5-5861-03	Medium Adult	3	195 mm	11 mm	1
5-5861-04	Large Adult	4	205 mm	16 mm	1
5-0234-21	Replacement Lamp	0, 1	—	—	1
5-0234-14	Replacement Lamp	2, 3, 4	—	—	1