


# UNI-LIM™ AND PEDI-LIM™ ANESTHESIA CIRCUIT

LATEX FREE | SINGLE PATIENT USE | DISPOSABLE

The single-limb circuit is the optimal choice for your Co-Axial Circle System needs. Designed to function as a standard circle anesthesia breathing system with streamlined tube-within-a-tube and eliminating bulky connectors.

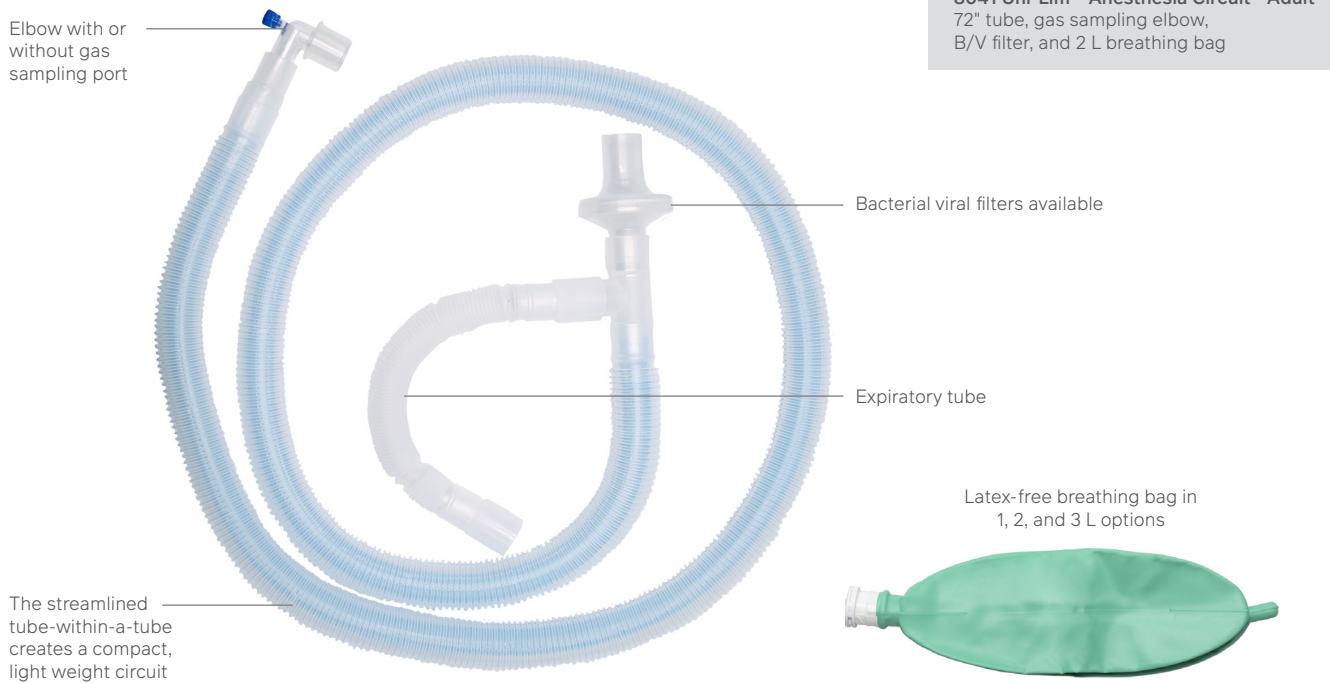


**8041 Uni-Lim™ Anesthesia Circuit - Adult**  
72" tube, gas sampling elbow, B/V filter, and 2 L breathing bag

ITEM	DESCRIPTION	SIZE	ELBOW	FILTER	HME/HMEF	BAG	GSL	MASK
9061	60" Uni-Lim	Adult	GSE	(2) B/V	—	2 L	10' x .050" ID M/M	Large Adult 900 Series
8041	72" Uni-Lim	Adult	GSE	B/V	—	2 L	—	—
9432	96" Uni-Lim	Adult	GSE	—	Maxi HMEF	—	—	—
6754	60" PEDI-Lim	Pediatric	GSE	(2) B/V	—	2 L	10' x .050" ID M/M	—

## UNI-LIM™ AND PEDI-LIM™ ANESTHESIA CIRCUITS

- Features a streamlined tube-within-a-tube design which eliminates twisting between the inspiratory and expiratory limbs
- Eliminates bulky connectors at the patient interface
- Provides natural heat and moisture conservation properties to that of a standard circle circuit
- Engineered with a secure-t inspiratory limb to prevent patient end disconnects – no floating internal tubes
- ISO-standard connections for secure, leak-resistant attachments

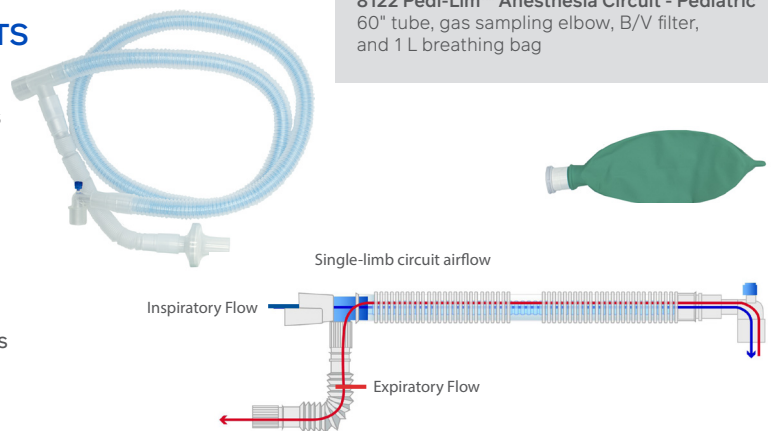

**8122 PEDI-Lim™ Anesthesia Circuit - Pediatric**  
60" tube, gas sampling elbow, B/V filter, and 1 L breathing bag

## CUSTOMIZE YOUR CIRCUIT

Choose from the following components to create a custom circuit to best meet your needs:




### BAGS

Three latex-free standard bag size options available (1, 2, and 3 L) and a 2 liter EZ gripper bag that is easy to grip with a gloved hand



### GAS SAMPLING LINES

Standard and swivel lines available

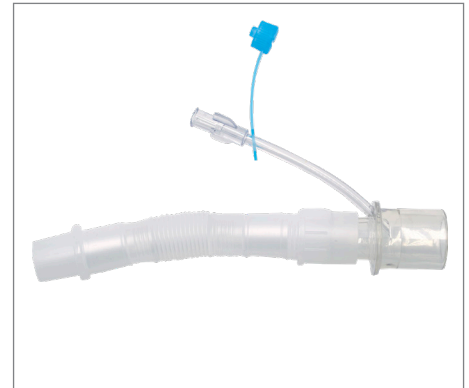

### FLEX EXTENDERS

5.2" Flex Tubing with and without capped gas sampling port



### FACE MASKS

Unique ergonomic mask shape that fits most facial configurations for improved sealing capabilities




### WYES AND CONNECTORS

Elbows, parallel and swivel wyes available with and without port



### FILTERS

HMEs, HMEFs, bacterial viral, and HEPA filters with luer gas sampling ports available. Sampling ports located at the top of the housing, just as if a common gas sampling elbow were being used. ISO-standard connections for secure, leak-resistant attachment to breathing systems.

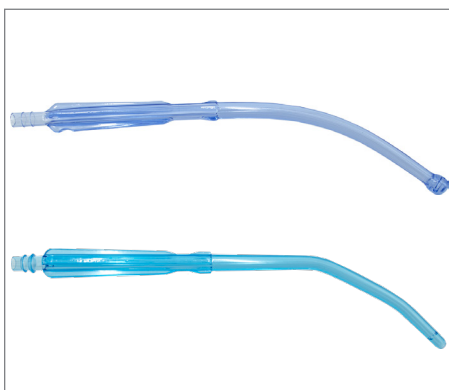

### YANKAUER

Sterile Yankauer with bulb tip and non-sterile Yankauer with regular tip available