


# T-PIECE NEONATAL PATIENT CIRCUIT

LATEX FREE | DEHP FREE | SINGLE PATIENT USE | DISPOSABLE

The T-piece resuscitator kit is a single patient use resuscitator circuit intended to provide respiratory support for neonatal patients. It is designed with a positive end-expiratory pressure (PEEP) valve to help maintain optimal pressure and made of corrugated tubing without any latex or added DEHP.

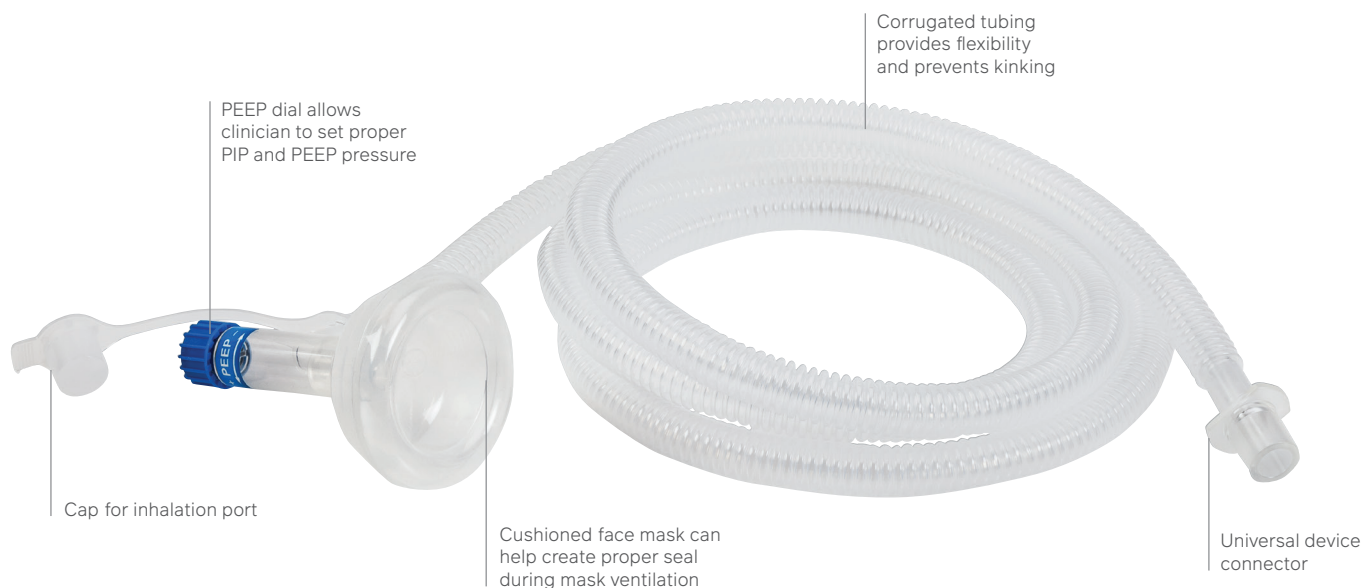

## T-PIECE RESUSCITATION KITS

- Universal device connector fits both 15 mm male and 10 mm female connections without the need for an additional adaptor
- Ventilation is provided by occluding PEEP dial hole intermittently at the desired respiratory rate. This allows for inspiration and exhalation to occur at the proper PIP and PEEP pressure
- Corrugated tubing is 6 feet long, enabling extended connection to the infant t-piece resuscitator
- Compatible with GE Giraffe Stand Alone Infant Resuscitation System
- Cushioned face mask available, size 0 mask has a 48 mm diameter, and a size 1 mask has a 60 mm diameter

ITEM	DESCRIPTION	MASK SIZE	PK
M1091335VS	T-Piece Neonatal Circuit	N/A	10
M1091316VS	T-Piece Neonatal Circuit Kit	Size 0	10
M1091365VS	T-Piece Neonatal Circuit Kit	Size 1	10