

# INCENTIVE SPIROMETERS


LATEX FREE | SINGLE PATIENT USE | DISPOSABLE

Deep breathing therapy is simplified with AirLife™ volumetric incentive spirometers. Their intuitive design prompts patients to correctly perform and monitor their own breathing exercises, even without direct supervision.

Dual-sided calibrations identifies achieved volumes up to 4,000 mL capacity. Clinicians and patients can easily read results from either side of the unit.

Easy to grip handle and helps ensure that patients use the spirometer in an upright position for proper therapy.

Adjustable indicator provides a volumetric goal for adequate lung exercise.

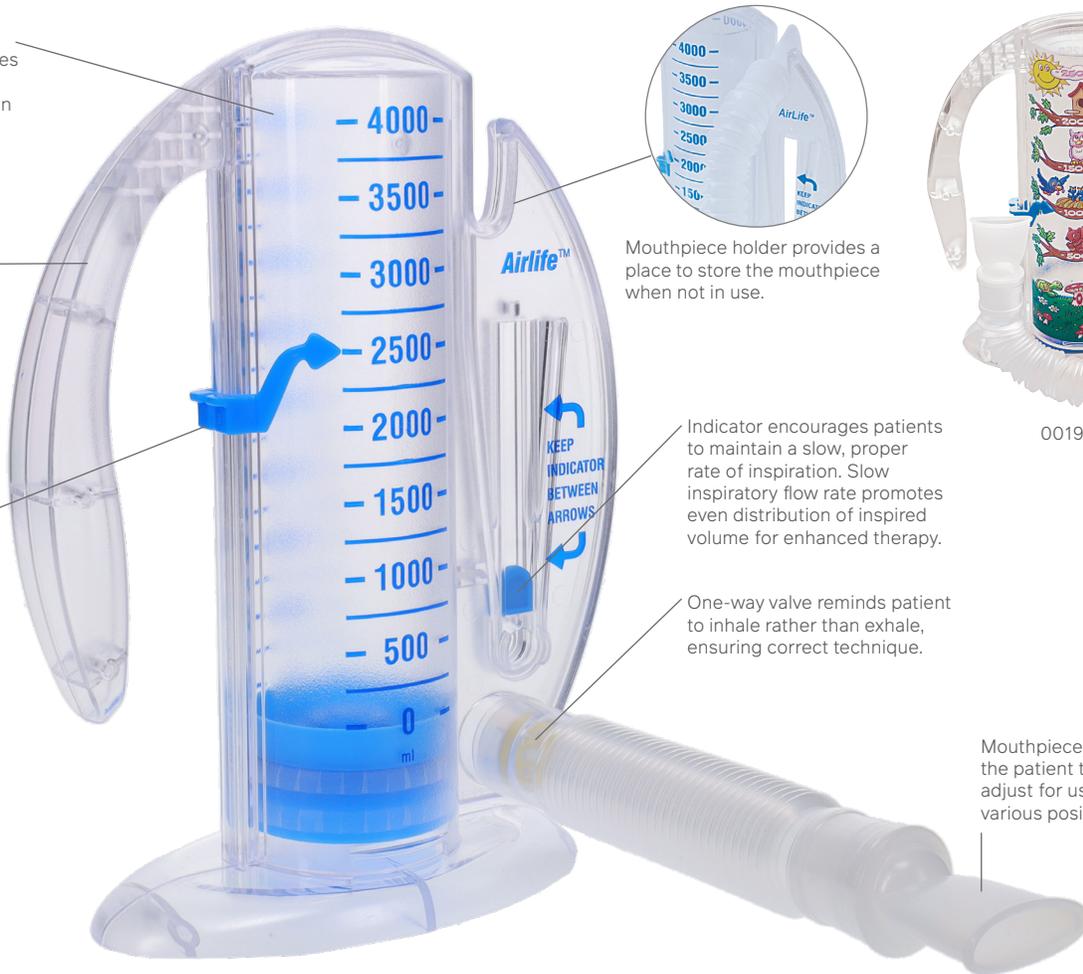

Mouthpiece holder provides a place to store the mouthpiece when not in use.



001905A

Indicator encourages patients to maintain a slow, proper rate of inspiration. Slow inspiratory flow rate promotes even distribution of inspired volume for enhanced therapy.

One-way valve reminds patient to inhale rather than exhale, ensuring correct technique.

Mouthpiece allows the patient to easily adjust for use in various positions.

## FEATURES

- The convenient handle can be placed on the bed rail for safe storage and easy accessibility.
- Durable plastic design resists accidental breakage for improved safety and hospital cost

ITEM	CAPACITY	SIZE	VALVE	PK
001901A	4,000 mL	Adult	One-Way Valve	12
001902A	4,000 mL	Adult	—	12
001903A	2,500 mL	Adult	One-Way Valve	12
001904A	2,500 mL	Adult	—	12
001905A	2,500 mL	Pediatric	One-Way Valve	12