

INFUSEIT™ PRESSURE INFUSOR

LATEX FREE | DEHP FREE | MR CONDITIONAL

InfuseIT™ pressure infusors are designed for ease of use. Faster inflation, individual hooks for solution and IV poles, and highly visible pressure gauge facilitate fast setup and easy maintenance. The enclosed valve eliminates the risk of fluids being trapped and hidden and the MRI conditional spring makes it safe to use within an MRI environment. Available in disposable and extended use reusable versions.



Palm-fitted, efficient hand bulb makes inflation easy.




Stopcock Valve


Highly visible pressure gauge provides easy to read 360° window which can be read from a distance



0-150 mm Hg

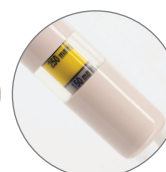

150-250 mm Hg

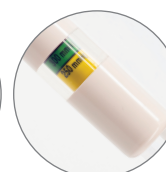

250-300 mm Hg



300+ mm Hg



Clear back panel provides full visibility of solution bag to monitor fluid level and easily scan patient barcode.



IV fluid bag hook is durable and designed to fit all standard bags snugly

Reusable bags are navy color and disposable bags are white.

FEATURES

- Ergonomically designed hand bulb for fast inflation
- Pressure relief valve prevents over inflation (375 mm Hg pressure relief)
- Enclosed valve helps eliminate the risk of fluids being trapped and hidden
- Pressure gauge utilizes a non-ferrous bronze spring making it MRI conditional
- Reusable InfuseIT models constructed of nylon for extended use
- Incredibly reliable - 100% tested

INFUSEIT™ PRESSURE INFUSOR

ITEM	DESCRIPTION	VOLUME	PK
ZIT-500-10	Disposable Pressure Infusor	500 mL	10
ZIT-1000-10	Disposable Pressure Infusor	1000 mL	10
ZIT-3000-10	Disposable Pressure Infusor	3000 mL	10
ZIT-520-5	Reusable Pressure Infusor	500 mL	5
ZIT-1020-5	Reusable Pressure Infusor	1000 mL	5